



Playback Switch

USER QUICK GUIDE

Introduction

The Playback Switch is an assistive switch that also plays a prerecorded message when the button is pressed. It has a similar functionality to the BIGmack and LITTLEmack buttons.

Features

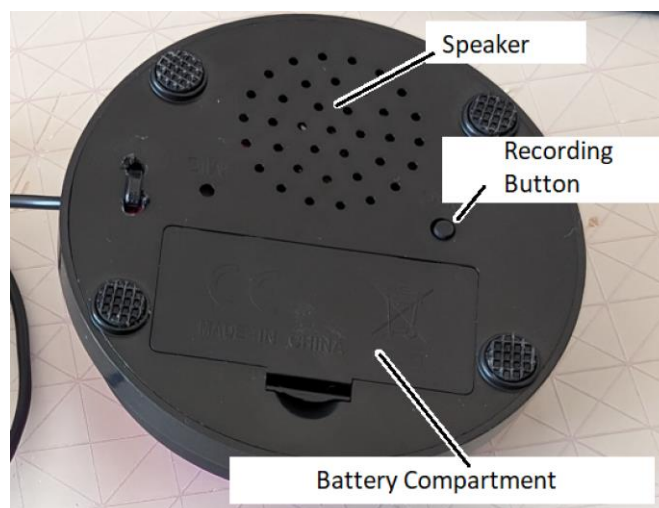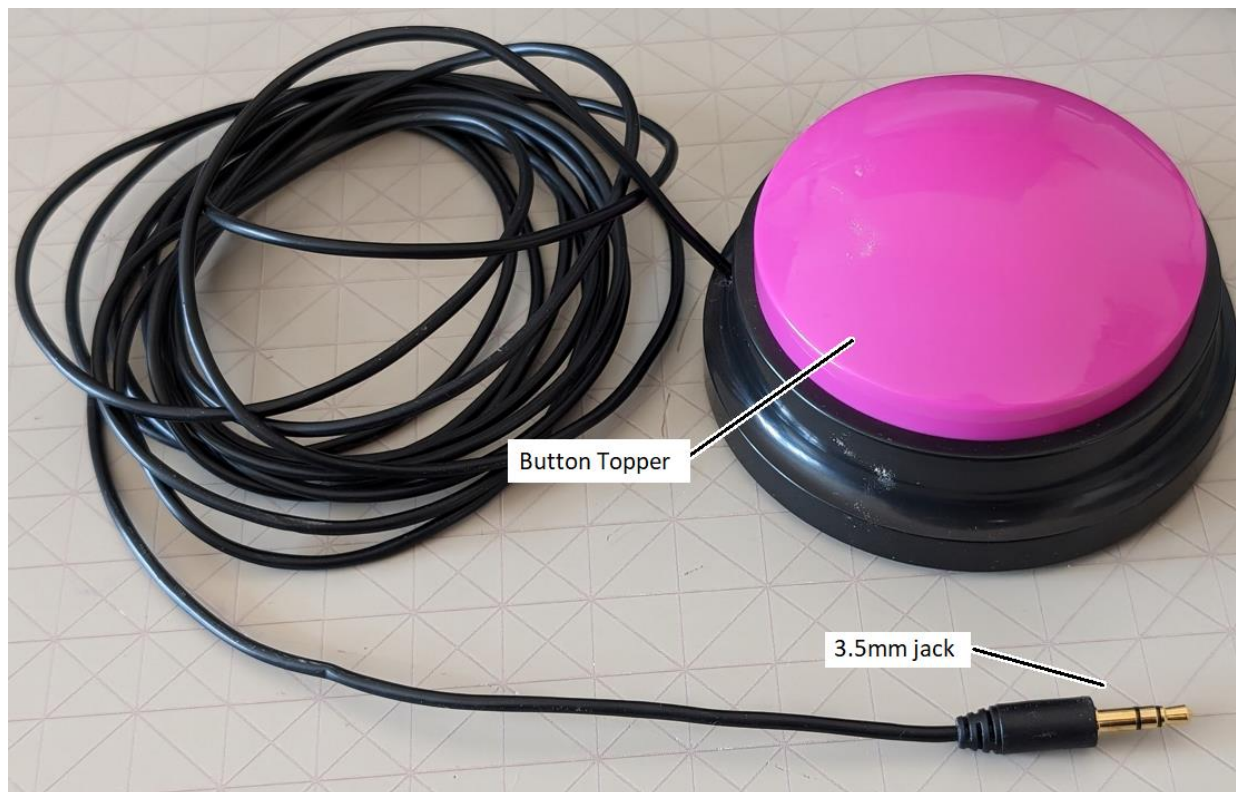

Playback Switch

USER QUICK GUIDE



Usage

To record a message, flip over the button and hold the small black button on the bottom. The button will beep and then record audio for 30 seconds or until the button is released, whichever comes first.

The 3.5mm cable can be connected to an accessible device, and then when the button is pressed it will activate the connected device and play back the recorded message.

Compatibility

This device is compatible with any accessible device that works with 3.5mm switches.

Specifications

Diameter: 90mm

Height: 45mm

Weight: 322 grams

Cleaning

Device can be wiped clean with a damp cloth.